



RECYCLING



POLYMER

PRODUCTSHEET

Intex Ortho




back

Properties




ADHESIVE

Both sides / Polyester
in compound

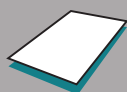

DELIVERY

1.00 x 1.50 m / 0.75 x 1.00 m

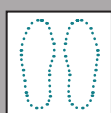

CUTTING

Length direction

A thermoplastic with an "anti-slip effect".

Thermoplastic and firm material covered with a tear proof polyester fabric on both sides. The edges turn fluid during processing and therefore end up smooth


Usage: for grip layer insoles.

Material Data

Types	Thickness	Size
Intex Ortho Anthracite	0,4 mm	1.00 x 1.50 m
Intex Ortho Anthracite	0,4 mm	0,75 x 1,00 m

Processing

Interface temperature
80°C-90°C

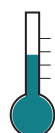

Activating temperature
130°C-140°C

5-8 seconds

Colors



Anthracite

For medical device Class 1
According (EU) MDR 2017/745



Reinforcements

Manufacturer's warranty & storing: When working with multiple materials during processing, chemical reactions and incompatibilities between components should always be tested beforehand. We are not liable for damage caused by failure to follow our instructions or responses beyond our influence and control. Under normal conditions, storage is unlimited, but avoid storage in direct sun. The goods are guaranteed for 6 months from the date of delivery.

www.frecoma.nl

Version 03/21